

Serene Paramount Visco



HCPCS: E2607/08 Skin Protection And Positioning Cushion

THE PARAMOUNT VISCO was designed for users that require moderate to high positioning support and provides excellent pressure relief. The cushion has aggressive leg troughs and abductor contouring that provides pelvic stability and leg positioning. A cushion made from medium density, high resilient, polyurethane foam with high density side wedges that provide rear side support. The paramount Visco also has 1" Visco insert which provide additional support to the IT area.

PRODUCT HIGHLIGHTS:

- ✓ Standard waterproof cover
- ✓ Radiused front edge for foot propelling
- ✓ Light weight cushion
- ✓ Nonslip bottom on washable , incontinent cover
- ✓ Promotes high postural support
- ✓ 4 way stretch, water resistant inner liner
- ✓ Medium density, high resilient, polyurethane foam
- ✓ High density polyurethane foam side wedges
- ✓ 1" high density visco foam insert


Product Weight: 3 Lbs / 1.36 Kg (18"W x20"L)




Non slip Bottom on washable, incontinent cover



Polyurethane foam, anatomically shaped

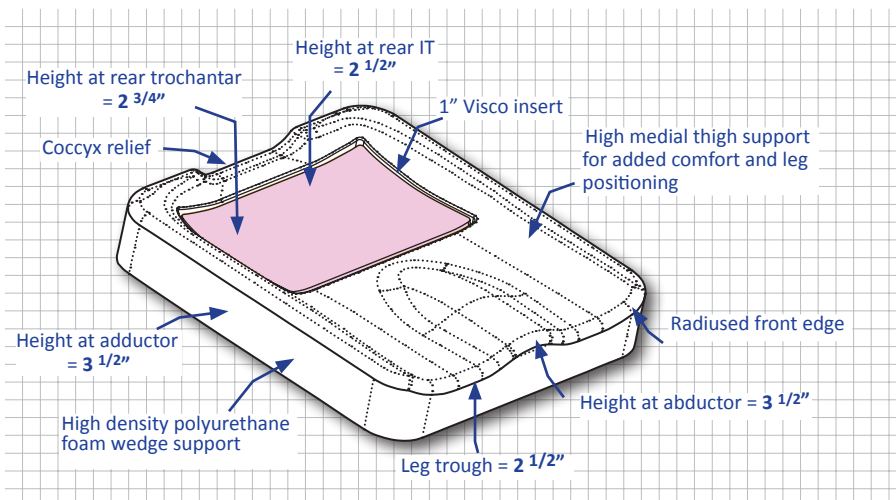

Pressure Reading Scale in mmHG

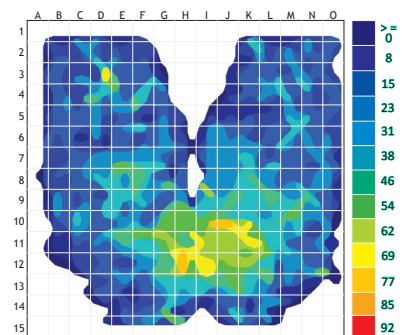

User Weight 250 Lbs / 113 Kg
 Duration 1 hr
 Average Pressure 31 mmHg

SERENE CUSHION STYLE: PARAMOUNT VISCO
Stock Sizes – Same Day Shipment

Part #	Size W x L	Part #	Size W x L
PAV1616	16" x 16"	PAV1818	18" x 18"
PAV1618	16" x 18"	PAV2016	20" x 16"
PAV1816	18" x 16"	PAV2018	20" x 18"

SERENE CUSHION STYLE: PARAMOUNT VISCO
Non-Stock Sizes – Shipment within 3 business days

Part #	Size W x L	Part #	Size W x L	Part #	Size W x L
PAV1620	16" x 20"	PAV2216	22" x 16"	PAV2416	24" x 16"
PAV1820	18" x 20"	PAV2218	22" x 18"	PAV2418	24" x 18"
PAV2020	20" x 20"	PAV2220	22" x 20"	PAV2420	24" x 20"