



# Supreme Gel Cushion

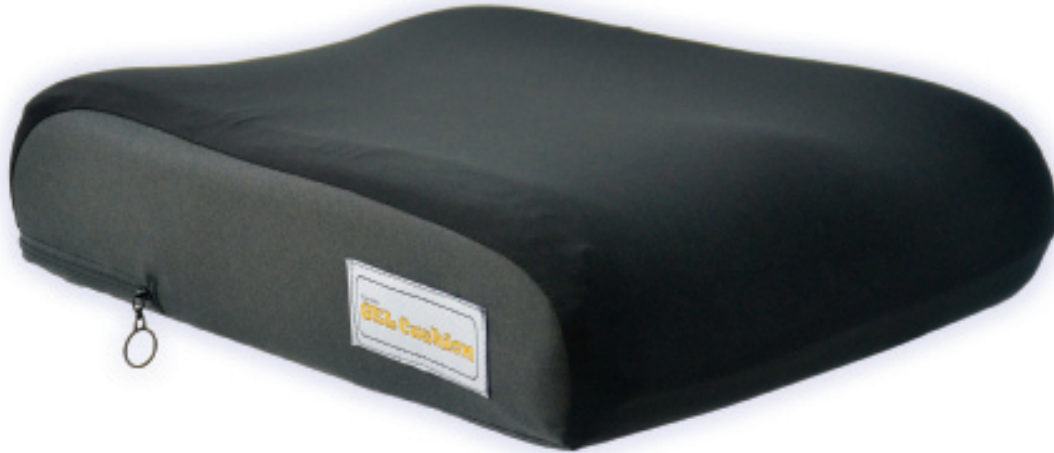

The Prism Supreme Gel™ Cushion was designed for users requiring mild to moderate positioning support and provides excellent pressure distribution relief. The gel overlay reduces shear and friction that contributes to skin breakdown, while the anatomically designed firm polyurethane base provides superior pelvic stability and leg positioning. A central core of pressure relieving visco-elastic insert under the IT's, allows for excellent immersion for bony prominences.

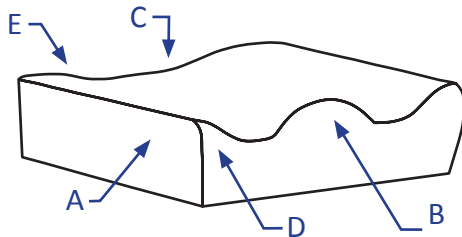

#### CUSHION THICKNESS

- A- Height at adductor = 5"
- B- Height at abductor = 5"
- C- Height at rear IT = 3 1/4"
- D- Leg Trough = 3 1/4"
- E- Height at rear trochanter = 4 1/4"

#### Prism Series Cushion STYLE: Supreme Gel

Stock Sizes – Same Day Shipment			
Part Numb	Size W x L	Part Number	Size W x L
SGC1616	16" x 16"	SGC1818	18" x 18"
SGC1618	16" x 18"	SGC1820	18" x 20"
SGC1816	18" x 16"		

Non-Stock Sizes – Shipment within 3 business days			
Part Number	Size W x L	Part Number	Size W x L
SGC1516	15" x 16"	SGC1916	19" x 16"
SGC1518	15" x 18"	SGC1918	19" x 18"
SGC1620	16" x 20"	SGC1920	19" x 20"
SGC1716	17" x 16"	SGC2016	20" x 16"
SGC1718	17" x 18"	SGC2018	20" x 18"
SGC1720	17" x 20"	SGC2020	20" x 20"

Custom Sizes - Shipment within 7 business days	
Part Number	Size W x L
SGCUS	Custom

**HCPCS E2607 / E2608**


# Supreme Gel Cushion



Cushion Layers and Inner Cover



## Product Highlights:

- ✓ Gel overlay provides even distribution of weight across the cushion
- ✓ Waterfall front edge for foot propelling comfort
- ✓ Cane cut in rear allowing correct placement in wheelchair and proper seat depth
- ✓ 3" of visco foam under IT's for superior immersion
- ✓ Supreme pressure relief
- ✓ Standard 4-way stretch, water resistant inner liner
- ✓ Standard 4-way stretch, breathable outer cover
- ✓ **Product Weight:** 12.5 lbs (18"W x 18"L)  
**User Weight:** 250 lbs  
**Warranty:** Two (2) years limited on cushion and 90 days on the cover.

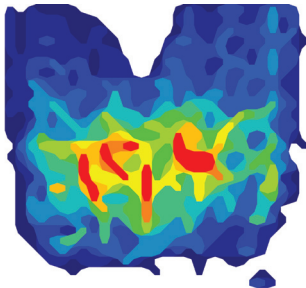

Pressure Map	Supreme Gel
Weight	220 Lbs/99 Kg
Duration	1 hr
Average Pressure	35 mmHg

Pressure Reading Scale in mmHG

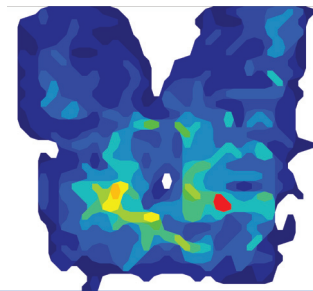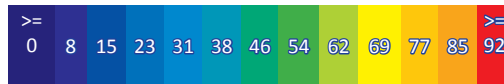

Pressure Map	Supreme Gel
Weight	150 Lbs/68 Kg
Duration	1 hr
Average Pressure	24 mmHg

